



## DRG<sup>®</sup> Adrenaline ELISA (EIA-4306)



Revised 12 Apr. 2011 rm (Vers. 8.1)

**RUO** in the USA

*Please use only the valid version of the package insert provided with the kit.*

*This kit is intended for Research Use Only.*

*Not for use in diagnostic procedures.*

### PRINCIPLE OF THE TEST

Enzyme Immunoassay for determination of Adrenaline (Epinephrine) in plasma and urine.

Adrenaline (epinephrine) is extracted by using a cis-diol-specific affinity gel, acylated and then derivatized enzymatically. The competitive ELISA kit uses the microtiter plate format. The antigen is bound to the solid phase of the microtiter plate. The derivatized standards, controls and samples and the solid phase bound analytes compete for a fixed number of antiserum binding sites. After the system is in equilibrium, free antigen and free antigen-antiserum complexes are removed by washing. The antibody bound to the solid phase is detected by an anti-rabbit IgG-peroxidase conjugate using TMB as a substrate. The reaction is monitored at 450 nm.

Quantification of unknown samples is achieved by comparing their absorbance with a reference curve prepared with known standard concentrations.

### 1. ADVICE ON HANDLING THE TEST

#### Reliability of the test results

In order to assure a reliable evaluation of the test results it must be conducted according to the instructions included and in accordance with current rules and guidelines (GLP, RILIBÄK, etc.). Special attention must be paid to control checks for precision and correctness during the test; the results of these control checks have to be within the norm range. In case of significant discrepancies between the pre-set assay characteristics of this test and the actual results please contact the manufacturer of the test kit for further instructions.

It is recommended that each laboratory establishes its own reference intervals. The values reported in this test instruction are only indicative.

#### Complaints

In case of complaints please submit to the manufacturer a written report containing all data as to how the test was conducted, the results received and a copy of the original test printout. Please contact the manufacturer to obtain a reclamation form and return it completely filled in to the manufacturer.

#### Warranty

This test kit was produced according to the latest developments in technology and subjected to stringent internal and external quality control checks. Any alteration of the test kit or the test procedure as well as the usage of reagents from different charges may have a negative influence on the test results and are therefore not covered by warranty. The manufacturer is not liable for damages incurred in transit.

#### Disposal

Residual substances and/or all remaining chemicals, reagents and ready for use solutions, are special refuse. The disposal is subject to the laws and regulations of the federation and the countries. About the removal of special refuse the



Revised 12 Apr. 2011 rm (Vers. 8.1)

**RUO** in the USA

responsible authorities or refuse disposal enterprises inform. The disposal of the kit must be made according to the national official regulations. Legal basis for the disposal of special refuse is the cycle economic- and waste law. The appropriate safety data sheets of the individual products are available upon request. The safety data sheets correspond to the standard: ISO 11014-1.

**Interference**

Do not mix reagents and solutions from different lots. Consider different transport and storage conditions. Inappropriate handling of test samples or deviations from the test regulation can the results affect. Use no kit components beyond the expiration date. Avoid microbiological contamination of the reagents and the washing water. Consider incubation periods and wash references.

**Precautions**

Observe the incubation periods and washing instructions. Never pipette by mouth and avoid contact of reagents and specimens with skin. No smoking, eating or drinking in areas where samples or kit test tubes are handled. When working with kit components or samples, always wear protective gloves and wash your hand thoroughly as soon as you have finished the work. Avoid spraying of any kind. Avoid any skin contact with reagents. Use protective clothing and disposable gloves. All steps have to be performed according to the protocol. Optimal test results are only obtained when using calibrated pipettes. Sodium azide could react with lead and copper tubes and may form highly explosive metal azide. When clearing up, rinse thoroughly with large volumes of water to prevent such formation.

All reagents of this testkit which contain human or animal serum or plasma have been tested and confirmed negative for HIV I/II, HbsAg and HCV by FDA approved procedures.

All reagents, however, should be treated as potential biohazards in use and for disposal.

**STORAGE AND STABILITY**

Store the reagents at 2 - 8 °C until expiration date.

Do not use components beyond the expiration date indicated on the kit labels.

Do not mix various lots of any kit component within an individual assay.

DRG® Adrenaline ELISA (EIA-4306)



Revised 12 Apr. 2011 rm (Vers. 8.1)

**RUO** in the USA

CONTENTS OF THE KIT

<b>FOILS</b>	<b>Adhesive Foil</b>	1 x 4	ready for use
<b>WASH-CONC</b> 50x	<b>Wash Buffer Concentrate</b>	1 x 20 mL	Concentrate. Dilute content with dist. water to a final volume of 1000 mL
<b>CONJUGATE</b>	<b>Enzyme Conjugate</b>	1 x 12 mL	ready for use, anti-rabbit IgG conjugated with peroxidase
<b>SUBSTRATE</b>	<b>Substrate</b>	1 x 12 mL	ready for use, containing a solution of tetramethylbenzidine (TMB)
<b>STOP-SOLN</b>	<b>Stop Solution</b>	1 x 12 mL	ready for use, containing 0.25 M H <sub>2</sub> SO <sub>4</sub>
<b>ADR</b> <b>ADR</b> <b>MN</b>	<b>Adrenaline-Metanephrine Microtiter Strips</b>	1 x 96 wells	12 strips, 8 wells each, break apart, pre-coated, blue coloured
<b>ADR-AS</b>	<b>Adrenaline Antiserum</b>	1 x 6 mL	from rabbit, ready for use, blue coloured, blue screw cap
<b>ADJUST-BUFF</b>	<b>Adjustment Buffer</b>	1 x 4 mL	ready for use
<b>STANDARD</b> <b>A</b>	<b>Standard A</b>	1 x 4 mL	ready for use
<b>STANDARD</b> <b>B</b>	<b>Standard B</b>	1 x 4 mL	ready for use
<b>STANDARD</b> <b>C</b>	<b>Standard C</b>	1 x 4 mL	ready for use
<b>STANDARD</b> <b>D</b>	<b>Standard D</b>	1 x 4 mL	ready for use
<b>STANDARD</b> <b>E</b>	<b>Standard E</b>	1 x 4 mL	ready for use
<b>STANDARD</b> <b>F</b>	<b>Standard F</b>	1 x 4 mL	ready for use
<b>ACYL-BUFF</b>	<b>Acylation Buffer</b>	1 x 20 mL	ready for use
<b>ACYL-REAG</b>	<b>Acylation Reagent</b>	1 x 3 ml	ready for use
<b>ASSAY-BUFF</b>	<b>Assay Buffer</b>	1 x 6 mL	ready for use, contains 1 M HCl
<b>COENZYME</b>	<b>Coenzyme</b>	1 x 2 mL	ready for use, S-adenosyl-L-methionine
<b>ENZYME</b>	<b>Enzyme</b>	2 x 1 mL	lyophilized, contains the enzyme COMT
<b>EXTRACT-BUFF</b>	<b>Extraction Buffer</b>	1 x 6 mL	ready for use
<b>EXTRACT-PLATE</b> 48	<b>Extraction Plate</b>	2 x 48 wells	coated with boronate affinity gel
<b>HCL</b>	<b>Hydrochloric Acid</b>	1 x 20 mL	ready for use, yellow coloured, contains 0.025 M HCl
<b>CONTROL</b> 1	<b>Control 1</b>	1 x 4 mL	ready for use
<b>CONTROL</b> 2	<b>Control 2</b>	1 x 4 mL	ready for use



Revised 12 Apr. 2011 rm (Vers. 8.1)

**RUO** in the USA**Additional materials and equipment required but not provided in the kit**

- Calibrated variable precision micropipettes (e.g. 10-100 µL / 100-1,000µL)
- Microtiter plate washing device, Absorbent material (paper towel)
- ELISA reader capable of reading absorbance at 450 nm and 620 or 650 nm
- Shaker (shaking amplitude 3mm; approx. 600 rpm)
- Distilled water
- Vortex mixer

**SAMPLE COLLECTION AND STORAGE****Plasma**

EDTA-Plasma should be used. Do not use haemolytic or lipemic samples.  
Storage: up to 6 hours at 2 - 8°C, for longer period (up to 6 month) at - 20°C.  
Repeated freezing and thawing should be avoided.

**Urine**

Spontaneous or 24-hours urine, collected in a bottle containing 10-15 mL of 6 M HCl, should be used.  
Storage: for longer periods (up to 6 months) at -20°C.  
Repeated freezing and thawing should be avoided.  
Avoid exposure to direct sunlight.

**TEST PROCEDURE**

Allow all reagents to reach room temperature and mix thoroughly by gentle inversion before use.  
Duplicate determinations are recommended.

**Preparation of reagents****Wash buffer**

Dilute the 20 mL Wash buffer Concentrate with distilled water to a final volume of 1000 mL.  
Storage: up to 6 months 4–8°C

**Enzyme solution**

Reconstitute the content of the vial labelled 'Enzyme' with 1 mL distilled water and mix thoroughly.  
Add 0.3 mL of Coenzyme followed by 0.7 mL of Adjustment Buffer.  
The total volume of the enzyme solution is 2.0 mL.

**■** *The enzyme solution has to be prepared freshly prior to the assay (not longer than 10 - 15 minutes in advance).  
Discard after use!*



Revised 12 Apr. 2011 rm (Vers. 8.1)

RUO
-----

 in the USA

### Sample preparation, extraction and acylation

1. Pipette **10 µL** of **standards, controls, urine samples** and **300 µL** of **plasma samples** into the respective wells of the **Extraction Plate**.
2. Add **250 µL** of **distilled water** to the wells with **standards, controls** and **urine samples**.
3. Pipette **50 µL** of **Assay Buffer** into all wells
4. Pipette **50 µL** of **Extraction Buffer** into all wells
5. Cover plate with adhesive foil and incubate **30 min** at **RT** (20-25°C) on a shaker (approx. 600 rpm).
6. Remove the foil. Empty plate and blot dry by tapping the inverted plate on absorbent material.
7. Pipette **1 mL** of **Wash Buffer** into all wells. Incubate the plate for **5 min** at **RT** (20-25°C) on a shaker (approx. 600 rpm). Empty plate and blot dry by tapping the inverted plate on absorbent material.
8. Pipette another **1 mL** of **Wash Buffer** into all wells. Incubate the plate for **5 min** at **RT** (20-25°C) on a shaker (approx. 600 rpm). Empty plate and blot dry by tapping the inverted plate on absorbent material.
9. Pipette **150 µL** of **Acylation Buffer** into all wells.
10. Pipette **25 µL** of **Acylation Reagent** into all wells.
11. Incubate **15 min** at **RT** (20-25°C) on a shaker (approx. 600 rpm).
12. Empty plate and blot dry by tapping the inverted plate on absorbent material.
13. Pipette **1 mL** of **Wash Buffer** into all wells. Incubate the plate for **10 min** at **RT** (20-25°C) on a shaker (approx. 600 rpm). Empty plate and blot dry by tapping the inverted plate on absorbent material.
14. Pipette **150 µL** of **Hydrochloric Acid** into all wells.
15. Cover plate with adhesive foil. Incubate **10 min** at **RT** (20-25°C) on a **shaker** (approx. 600 rpm). Remove the foil and discard.

#### **Do not decant the supernatant thereafter!**

The following volumes of the supernatant are needed for the subsequent ELISA:

<b>Adrenaline</b>	<b>100 µL</b>
-------------------	---------------

### Adrenaline ELISA

1. Pipette **25 µl** of the **Enzyme Solution** (refer to 6.1) into all wells of the **Adrenaline Microtiter Strips**.
2. Pipette **100 µL** of the extracted **standards, controls and samples** into the appropriate wells.
3. Incubate for **30 min** at **RT** (20-25°C) on a shaker (approx. 600 rpm).
4. Pipette **50 µL** of the respective **Adrenaline Antiserum** into all wells and cover plate with **Adhesive Foil**.
5. Incubate for **2 hours** at **RT** (20-25°C) on a shaker (approx. 600 rpm).
6. Remove the foil. Discard or aspirate the content of the wells and **wash** each well **3 times** thoroughly with **300 µL Wash Buffer**. Blot dry by tapping the inverted plate on absorbent material.



Revised 12 Apr. 2011 rm (Vers. 8.1)

**RUO** in the USA

7. Pipette **100 µL** of the **Enzyme Conjugate** into all wells.
8. Incubate for **30 min** at **RT** (20-25°C) on a shaker (approx. 600 rpm).
9. Discard or aspirate the content of the wells and **wash** each well **3 times** thoroughly with **300 µL Wash Buffer**. Blot dry by tapping the inverted plate on absorbent material.
10. Pipette **100 µL** of the **Substrate** into all wells and incubate for **25 ±5 min** at **RT** (20-25°C) on a shaker (approx. 600 rpm). **■ Avoid exposure to direct sun light!**
11. Add **100 µL** of the **Stop Solution** to each well and shake the microtiter plate to ensure a homogeneous distribution of the solution
12. **Read** the absorbance of the solution in the wells within 10 minutes, using a microplate reader set to **450 nm** and a reference wavelength between 620 nm and 650 nm.

**CALCULATION OF RESULTS**

	Concentration of the standards					
Standard	A	B	C	D	E	F
Adrenaline (ng/mL)	0	1	4	15	50	200
Adrenaline (nmol/L)	0	5.5	22	82	273	1 092
Conversion:	Adrenaline (ng/mL) x 5.46 = Adrenaline (nmol/L)					

The calibration curve is obtained by plotting the absorbance readings (calculate the mean absorbance) of the standards (linear, y-axis) against the corresponding standard concentrations (logarithmic, x-axis). Use a non-linear regression for curve fitting (e.g. spline, 4- parameter, akima).

**Urine samples and controls:**

The concentrations of the **urine samples** and the **Controls 1 & 2** can be read directly from the standard curve. Calculate the 24 h excretion for each urine sample: **µg/24h = µg/L x L/24h**

**Plasma samples:**

The read concentrations of the **plasma samples** have to be **divided by 30**.

**Quality control**

It is recommended to use control samples according to state and federal regulations. Use controls at both normal and pathological levels. The kit, or other commercially available, controls should fall within established confidence limits. The confidence limits of the kit controls are printed on the QC-Report.

**Calibration**

The binding of the antisera and the enzyme conjugates and the activity of the enzyme used are temperature dependent, and the extinction values may vary if a thermostat is not used. The higher the temperature, the higher the extinction values will be. The extinction values also depend on the incubation times. The optimal temperature during the Enzyme Immunoassay is between 20-25°C.

**■ In case of overflow, read the absorbance of the solution in the wells within 10 minutes, using a microplate reader set to 405 nm**




Revised 12 Apr. 2011 rm (Vers. 8.1)

**RUO** in the USA

### Typical calibration curve

*Example. Do not use for calculation!*

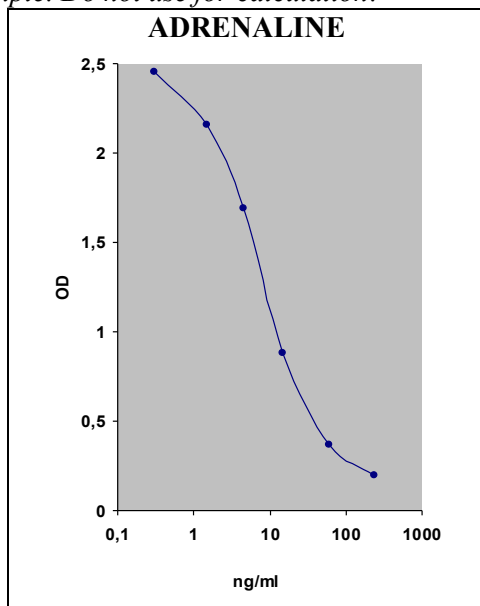

Version 2011-03-31~rm